

DemoThermometryDICOM

November 13, 2021

1 Demo of Proteus MR Thermometry

Redistribution and use in source and binary forms, with or without modification, are permitted provided that the following conditions are met:

- * Redistributions of source code must retain the above copyright notice, this list of conditions and the following disclaimer.
- * Redistributions in binary form must reproduce the above copyright notice, this list of conditions and the following disclaimer in the documentation and/or other materials provided with the distribution.
- * Neither the name of the University of Calgary nor the names of its contributors may be used to endorse or promote products derived from this software without specific prior written permission.

THIS SOFTWARE IS PROVIDED BY THE COPYRIGHT HOLDERS AND CONTRIBUTORS "AS IS" AND ANY EXPRESS OR IMPLIED WARRANTIES, INCLUDING, BUT NOT LIMITED TO, THE IMPLIED WARRANTIES OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE ARE DISCLAIMED. IN NO EVENT SHALL University of Calgary, Samuel Pichardo or any of the contributors be liable for any damages, INCLUDING DIRECT, INDIRECT, INCIDENTAL, SPECIAL, EXEMPLARY, OR CONSEQUENTIAL DAMAGES (INCLUDING, BUT NOT LIMITED TO, PROCUREMENT OF SUBSTITUTE GOODS OR SERVICES; LOSS OF USE, DATA, OR PROFITS; OR BUSINESS INTERRUPTION) HOWEVER CAUSED AND ON ANY THEORY OF LIABILITY, WHETHER IN CONTRACT, STRICT LIABILITY, OR TORT (INCLUDING NEGLIGENCE OR OTHERWISE) ARISING IN ANY WAY OUT OF THE USE OF THIS SOFTWARE, EVEN IF ADVISED OF THE POSSIBILITY OF SUCH DAMAGE.

This notebook illustrates the basic operation to use Proteus MR Thermometry library

Be aware some of the underlying structure for the processing is aligned how the main Proteus GUI application organizes the data.

```
[73]: %matplotlib inline
import matplotlib.pyplot as plt

import numpy as np
from Proteus.ThermometryLibrary import ThermometryLib
from Proteus.File_IO.H5pySimple import ReadFromH5py, SaveToH5py
import tables
import logging
from pprint import pprint
from Proteus.ThermometryLibrary.ThermometryLib import  LOGGER_NAME

from skimage import data, img_as_float
from skimage import exposure
import warnings
```

```

logger = logging.getLogger(LOGGER_NAME)

stderr_log_handler = logging.StreamHandler()
logger.addHandler(stderr_log_handler)

formatter = logging.Formatter('%(asctime)s - %(name)s - %(levelname)s -\n
→%(message)s')

logger.setLevel(logging.ERROR) #use logging.INFO or logging.DEBUG for detailed\n
→step information for thermometry

stderr_log_handler.setFormatter(formatter)

logger.info(': INITIALIZING MORPHEUS application, powered by Proteus')

```

1.1 Functions and classess required to prepare the MR processing

```

[74]: def CompareTwoOrderedLists(list1,list2,bPrintResults=False):
    '''
    Tool function to evaluate two MRI data collections are equivalent
    '''
    #please note NavigatorData will be empty, we kept for completeness purposes
    [IMAGES2,Navigator2]=list2
    [IMAGES,Navigator]=list1
    badimagecount = 0
    badimages = []
    badimageindex = []

    notallsame = False

    #Element wise comparison of two sets of results to ensure they match\n
    →eachother within tolerance
    for Stack in IMAGES:
        for Map in IMAGES[Stack]:
            if Map in ['TimeArrival','SelPointsROI', 'MaskROI'\n
            →,'TemperatureROIIMask']:
                continue
            for Number in range(len(IMAGES[Stack][Map])):
                for Slice in range(len(IMAGES[Stack][Map][Number])):

                    for Data in IMAGES[Stack][Map][Number][Slice]:
                        if type(IMAGES[Stack][Map][Number][Slice][Data]) is np.
                        →ndarray:
                            comparison=np.all(np.
                        →isclose(IMAGES[Stack][Map][Number][Slice][Data],

```

```

    ↪ IMAGES2[Stack][Map][Number][Slice][Data]))
        if comparison == False:
            notallsame=True
            if badimageindex.count(Number) == 0:
                badimagecount += 1
                badimageindex.append(Number)
            badimages.
    ↪ append((badimagecount,Stack,Map,Number,Slice,Data))
        elif type(IMAGES[Stack][Map][Number][Slice][Data]) is
    ↪ dict:
        for k in IMAGES[Stack][Map][Number][Slice][Data]:
            v1=IMAGES[Stack][Map][Number][Slice][Data][k]
            v2=IMAGES2[Stack][Map][Number][Slice][Data][k]
            if type(v1) is np.ndarray:
                comparison=np.all(np.isclose(v1,v2))
            else:
                comparison=v1==v2
            if comparison == False:
                notallsame=True
                if badimageindex.count(Number) == 0:
                    badimagecount += 1
                    badimageindex.append(Number)
                badimages.
    ↪ append((badimagecount,Stack,Map,Number,Slice,Data,k))
        else:
            comparison =
    ↪ (IMAGES[Stack][Map][Number][Slice][Data] ==
    ↪ IMAGES2[Stack][Map][Number][Slice][Data])
            if comparison == False:
                notallsame=True
                if badimageindex.count(Number) == 0:
                    badimagecount += 1
                    badimageindex.append(Number)
                badimages.
    ↪ append((badimagecount,Stack,Map,Number,Slice,Data))

    if bPrintResults:
        if notallsame == True:
            if len(badimages)>0:
                print ('The following images did not match within tolerance')
                for e in badimages:
                    print(e)
            else:

```

```

        print('*'*40+'\nDatasets were equivalent')

    return notallsame==False

def CreateSortedDataForProcessing(OBJ):
    '''
        The two main results to extract for processing are the images and navigator_
        ↪data dictionaries
        For thermometry processing , we only need to recover magnitude and phase_
        ↪data
    '''

    IMAGES=OBJ['IMAGES']
    NavigatorData=OBJ['ExtraData']['NavigatorData']
    IMAGES2 = {}
    ALL_ITEMS=[]
    for k in IMAGES:
        #this helps to initializes some empty data structures
        IMAGES2[k]={ 'MaskROI': [None], 'SelPointsROI': [None]}
        for k2 in { 'MaskROI': [None], 'SelPointsROI': [None]}:
            IMAGES2[k][k2] = IMAGES[k][k2]

    #we reorder the data to mimic how it comes when collecting from a real MRI_
    ↪scanner
    for SelKey in IMAGES:
        for StackMag,StackPhase in_
        ↪zip(IMAGES[SelKey]['Magnitude'],IMAGES[SelKey]['Phase']):
            for ImagMag,ImagPhase in zip(StackMag,StackPhase):
                ALL_ITEMS.append(ImagMag)
                ALL_ITEMS.append(ImagPhase)
    ALL_ITEMS.extend(NavigatorData)
    #the data is organized by time of arrival, to emulate how it works during_
    ↪MRI data collection
    ORDERED_ITEMS = sorted(ALL_ITEMS, key=lambda k: k['TimeStamp'])
    return IMAGES2,ORDERED_ITEMS

class InspectMPSData(object):
    '''
        Minimal class to open MPS files for the re processing
    '''
    def __init__(self,fname):
        print('fname',fname)
        self.fname=fname
        self.ATables=tables.open_file(fname,'r')
        A=self.ATables

```

```

NumberTreatments=A.root.Data.MRIONLINE._g_getnchildren()
print("Number of treatments ",NumberTreatments)

for treatment in A.root.Data.MRIONLINE._f_list_nodes():
    print('    '+treatment._v_name)

def GetDataTreatment(self,iddata):
    node=self.ATables.get_node("/Data/MRIONLINE/" +iddata)
    print(node)
    return ReadOnlineMRIData(node)

class FieldsForImaging:
    '''
    Class containing attributes defining the parameters controlling the
    →thermometry

    The values can be adjusted to test different parameter conditions
    '''
    def __init__(self):
        self.Alpha = 9.4e-09 #Thermometry temperature coefficient
        self.Beta = 3.0 # Beta Coefficient
        self.Gamma = 42580000.0 #Gyromagnetic ratio
        self.T_tolerance = 12.0 #SNR limit (*C)
        self.CorrectionOrder = 0 #Order of drift correction
        self.NumberOfAverageForReference = 4 #number of dynamics for averaging
        self.StartReference = 4 #dyn. index ref., thermometry is not
    →calculated in dynamics prior to this #
        self.TBaseLine = 37 #Baseline temperature
        self.CalculateLargeHistory = True #Calculate extra history
        self.UseUserDriftMask = True # use user-specified ROIs to select mask
    →for drift corrector
        self.ROIs = '1 C 4' # string defining ROI mask for monitoring, take a
    →look at \Proteus\Tools\parseROIString.py for details for use
        self.UserDriftROIs = '1 R 25 12 0 25' # string defining ROI mask for
    →drift corrector, take a look at \Proteus\Tools\parseROIString.py for details
    →for use
        #old mask settings for drift, better to UserDriftROIs instead
        self.CircleSizeFORSNRCoronal=45.0
        self.RectSizeFORSNRTransverse=110.0
        self.MaxSizeSNRRegion=200.0

        self.UseTCoupleForDrift = False #use this if have a setting using
    →thermocouples to minimize excessive drift correction

        self.NumberSlicesCoronal = 1 #Number of slices in coronal stack
        self.T_mask = 38.0 #Lower limit for temperature mask

```

```

#ECHO NAVIGATOR MOTION COMPENSATOR RELATED parameters
# just kept for completeness as they are now rarely used as we do not
→have anymore the echonavigator patch
self.UseMotionCompensation = False #Use Motion Compensation, keep this
→FALSE unless you have a dataset with echo navigator
self.TimeBeforeFilterNavigator = 10.0 #time before filtering (s)
self.OrderForPredictor = 5 #Order of predictor
self.DiscardedPointsInPredictor = 100 #Tail points to ignore
self.AmplitudeCriteriaForRestMotion = 25.0 # ampl. limit for
→motion-less detection (%)
self.TimeWindowForClassification = 11 #time window for class. (s)
self.TimeWindowForFiltering = 100 #time window for filter. (s)
self.NumberPointsInterpolateInitialLUT = 100 #Number of points for
→interpolation fir
self.NumberNavMessagesToWait = 0 #Number of Navigator messages to wait
→for
self.TimeWindowtoKeepInLUT = 175.0 #'Length of window (s) of entries
→to keep in LUT'
self.FrequencyCut = 0.8 #Frequency cutoff for butterworth filter (Hz)

#Empty Main object to preserve the structure required by thermometrylib
class MainObject: pass

class UnitTest:
    def __init__(self):
        #setting up supporting structures required to perform thermometry
        self.ImagingFields=FieldsForImaging()

        self.MainObject = MainObject()
        self.MainObject.TemporaryData = {}
        self.MainObject.TemporaryData['NavigatorDisplacement']=[]
        self.MainObject.TemporaryData['FilterForNavigator']=[]
        self.MainObject.NavigatorData=[]
        self.MainObject.ImagesKeyOrder=['Coronal', 'Sagittal', 'User1', 'User2']
        self.MainObject.IMAGES={}
        for k in self.MainObject.ImagesKeyOrder:
            self.MainObject.IMAGES[k]={'Magnitude': [], 'Phase': [], 'Temperature':
→[], 'Dose': [], 'MaskROI': [None], 'SelPointsROI': [None],
            'TemperatureROIMask': [None]}
            self.MainObject.TemporaryData[k]=[]
        self.POOL_SIZE=10000
        self.POOL_TIME_NAV=np.zeros(self.POOL_SIZE)
        self.POOL_DATA_NAV=np.zeros(self.POOL_SIZE)
        self.POOL_FILT_DATA_NAV=np.zeros(self.POOL_SIZE)

```

```

self.POOL_MOTIONLESS=np.ones(self.POOL_SIZE)*np.nan
self.POOL_INHALATION=np.ones(self.POOL_SIZE)*np.nan
self.POOL_EXHALATION=np.ones(self.POOL_SIZE)*np.nan
self.POOL_FILT_DATA_CLASS=np.zeros(self.POOL_SIZE)
self.POOL_DATA_INDEX=0

self.BackPointsToRefresh=200
self.TotalImages = 0
self.BottomIndexForFiltering=0
self.TProcessor={}
self.InBackground = False
self.cback_UpdateTemperatureProfile = lambda x: None
self.cback_UpdateNavigatorDisplacementProfile = lambda: None
self.cback_UpdateMRIImage = lambda x,y,z: None
self.cback_LockMutex = lambda x: None
self.cback_LockList = lambda x: None
self.IncreaseCounterImageProc = lambda: None
self.MaxSlicesPerDynamicProc = lambda: 1 #THIS IS ONLY VALID FOR DATA
↳ COLLECTIONS WITH 1 slice per dynamic
self.GetStackFromSliceNumberFunc = lambda x: (0,0)
self.NumberSlicesStackFunc = lambda x: 1 #THIS IS ONLY VALID FOR DATA
↳ COLLECTIONS WITH 1 slice per dynamic
self.ReleaseOnlyNavigatorProc = lambda: None

def ReturnElementsToInitializeprocessor(self):
    '''
    This function prepares a minimal
    '''
    MO=self.MainObject
    return [MO.IMAGES,
            MO.TemporaryData,
            MO.NavigatorData,
            MO.ImagesKeyOrder,
            self.IncreaseCounterImageProc,
            self.MaxSlicesPerDynamicProc,
            self.GetStackFromSliceNumberFunc,
            self.NumberSlicesStackFunc,
            self.ReleaseOnlyNavigatorProc]

def BatchProcessor(self, inputdata):
    '''
    This function reprocess all the magnitude and phase data to recreate a
    ↳ new data collection including thermometry data
    '''
    #add input data to parent

```

```

        [IMAGES2,self.MainObject.ORDERED_ITEMS] = _
    ↪CreateSortedDataForProcessing(inputdata)
    self.MainObject.MinTime = self.ep.GetReferenceTime()
    #process entries one by one
    for NewEntry in self.MainObject.ORDERED_ITEMS:

        if 'info' in NewEntry:
            self.ep.ProcessImage(NewEntry)
            self.TotalImages+=1

        else:
            self.ep.ProcessNavigator(NewEntry)
            self.TotalImages+=1

    return [self.MainObject.IMAGES,self.MainObject.NavigatorData]

def BatchProcessorFromList(self, ListInputdata):
    '''
    This function reprocess all the magnitude and phase data to recreate a
    ↪new data collection including thermometry data
    '''

    #add input data to parent
    self.MainObject.ORDERED_ITEMS=ListInputdata
    self.MainObject.MinTime = 0.0
    #process entries one by one
    for NewEntry in self.MainObject.ORDERED_ITEMS:

        if 'info' in NewEntry:
            self.ep.ProcessImage(NewEntry)
            self.TotalImages+=1

        else:
            self.ep.ProcessNavigator(NewEntry)
            self.TotalImages+=1

```

1.2 function to load DICOM

```

[117]: import sys
import glob
import os
import pydicom as dicom
def LoadDICOMGe(DCMDir='./',NumberSlicesCoronal =1,ManualTimeBetweenImages=5.0):
    '''
    Function to load an MR dataset from GE Scanner.

```



```

'''

AllFiles=glob.glob(DCMDir+os.sep+'*.dcm')
AllFiles.sort()
NumberFiles=len(AllFiles)

#if must be mutltiple of number of dynamics x 2 (real,imag)
if NumberFiles%3 !=0:
    raise ValueError('The number of images must be a multiple of 3')

#we'll scan and store all the images and verify how many stacks are in
→function of the different image orientation

MatPosOrientation=np.zeros((4,4))
MatPosOrientation[3,3]=1

AllImag=[]
am=None
ar=None
ai=None
pPos = None
nDynamic=1
PreSort=[]
for n in range(NumberFiles):
    fdcm = dicom.read_file(AllFiles[n])
    PreSort.append(fdcm)

warnings.warn('The DICOMS are missing TriggerTime in their Metadata, \n so
→there is no automatic way to recover the timing ')
#PreSort.sort(key=lambda fdcm: float(fdcm.TriggerTime))

#print(PreSort)

for fdcm in PreSort:
    if fdcm[0x0043, 0x102f].value==0:
        if am is not None:
            raise ValueError('There should not be preloaded magnitude')
        am = fdcm
    elif fdcm[0x0043, 0x102f].value==2:
        if ar is not None:
            raise ValueError('There should not be preloaded real')
        ar = fdcm
    elif fdcm[0x0043, 0x102f].value==3:
        if ai is not None:
            raise ValueError('There should not be preloaded imag')

```

```

        ai = fdcm
    else:
        raise ValueError('unhandled image type' +str(fdcm))

    if am is None or ai is None or ar is None:
        continue
    im = am
    if pPos is not None:
        if pPos==im.ImagePositionPatient:
            nDynamic+=1
    for m in range(2):
        entry={}
        entry['TimeStamp']=nDynamic*ManualTimeBetweenImages
        #float(im.TriggerTime)/1000.0
        S1={}
        if m==0:
            imdata=(im.pixel_array).astype(np.float32)
        else:
            cdata= (ar.pixel_array).astype(np.float32)+\
                (ai.pixel_array).astype(np.float32) *1j
            imdata=-np.angle(cdata)

        S1['VoxelSize']=np.zeros(3)
        S1['VoxelSize'][0:2]=np.array(im.PixelSpacing)/1e3
        S1['VoxelSize'][2]=float(im.SliceThickness)/1e3
        S1['DynamicLevel']=nDynamic
        S1['EchoTime']=float(im.EchoTime)/1e3
        S1['DynamicAcquisitionTime']=entry['TimeStamp']
        S1['ImageOrientationPatient']=np.array(im.ImageOrientationPatient)
        S1['ImagePositionPatient']=np.array(im.ImagePositionPatient)

        ImagePositionPatient=np.array(im.ImagePositionPatient)
        ImageOrientationPatient=np.array(im.ImageOrientationPatient)
        VoxelSize=np.array(im.PixelSpacing)
        MatPosOrientation=np.zeros((4,4))
        MatPosOrientation[3,3]=1
        MatPosOrientation[0:3,0]=ImageOrientationPatient[0:3]*VoxelSize[0]
        MatPosOrientation[0:3,1]=ImageOrientationPatient[3:]*VoxelSize[1]
        MatPosOrientation[0:3,3]=ImagePositionPatient

        CenterRow=im.Rows/2
        CenterCol=im.Columns/2
        IndCol=np.zeros((4,1))
        IndCol[0,0]=CenterRow
        IndCol[1,0]=CenterCol
        IndCol[2,0]=0
        IndCol[3,0]=1

```

```

CenterImagePosition=np.dot(MatPosOrientation,IndCol)

S1['OffcentreAnteriorPosterior']=CenterImagePosition[1,0]
S1['OffcentreFeetHead']=CenterImagePosition[2,0]
S1['OffcentreRightLeft']=CenterImagePosition[0,0]

S1['RescaleSlope']=1.0
S1['RescaleIntercept']=0.0
S1['ScaleSlope']=1.0
S1['ScaleIntercept']=0.0
S1['SlicePrepulseDelay']=0
S1['IsPhaseImage']=(m!=0)
NewEntry={'TimeStamp':entry['TimeStamp'],'info':S1,'data':imdata}
AllImag.append(NewEntry)

am=None
ar=None
ai=None
if pPos is None:
    pPos=im.ImagePositionPatient

#we recreated a pseudo-arrival by ordering the images by timestamp and by
→type (mag or phase)
SortedImag=sorted(AllImag, key=lambda d:
→(d['TimeStamp'],d['info']['IsPhaseImage']))

ListOfStacks=[]
FinalList=[]
for entry in SortedImag:
    S1=entry['info']['ImageOrientationPatient'].
→tolist()+entry['info']['ImagePositionPatient'].tolist()
    if S1 not in ListOfStacks:
        ListOfStacks.append(S1)
    for entry in SortedImag:
        S1=entry['info']['ImageOrientationPatient'].
→tolist()+entry['info']['ImagePositionPatient'].tolist()
        entry['info']['SliceNumber']=ListOfStacks.index(S1)
        FinalList.append(entry)
return FinalList

```

```

[151]: ListDataExample=LoadDICOMGe('MR_Thermometry_Examples/ExampleB/All/
→',ManualTimeBetweenImages=5.0)
print('Total number of images (both magnitude and phase)')
→',len(ListDataExample))
print('Basic Metatada')

```

```
pprint(ListDataExample[0]['info'])
```

Total number of images (both magnitude and phase) = 20

Basic Metadata

```
{'DynamicAcquisitionTime': 5.0,
  'DynamicLevel': 1,
  'EchoTime': 0.008400000000000001,
  'ImageOrientationPatient': array([ 9.99740e-01,  2.44158e-04, -2.27977e-02,
-2.27971e-02,
    1.09560e-02, -9.99680e-01]),
  'ImagePositionPatient': array([-145.303 , -46.4044, 139.583 ]),
  'IsPhaseImage': False,
  'OffcentreAnteriorPosterior': -44.8363061989888,
  'OffcentreFeetHead': -3.5704218572800075,
  'OffcentreRightLeft': -8.524741565439996,
  'RescaleIntercept': 0.0,
  'RescaleSlope': 1.0,
  'ScaleIntercept': 0.0,
  'ScaleSlope': 1.0,
  'SliceNumber': 0,
  'SlicePrepulseDelay': 0,
  'VoxelSize': array([0.0010938, 0.0010938, 0.01      ])}
```

/Users/spichardo/.edm/envs/MORPHEUS/lib/python3.6/site-packages/ipykernel_launcher.py:37: UserWarning: The DICOMS are missing TriggerTime in their Metadata, so there is no automatic way to recover the timing

We use the UnitTest class for demonstration to reprocess MRI data

```
[169]: ut = UnitTest() #Instantiate a parent class
ut.ep = ThermometryLib.EntryProcessing(*ut.
    ↪ReturnElementsToInitializeprocessor()) #Instantiate an entry processor
    ↪member on the parent class
ut.ep.ImagingFields = ut.ImagingFields #Instantiate a class full of image
    ↪processing parameters
ut.ep.ImagingFields.Beta=1.5
ut.ep.ImagingFields.ROIs='1 C 3 -13 -18'
ut.ep.ImagingFields.UserDriftROIs='-1 R 70 50 0 -10'
ut.ep.ImagingFields.UseUserDriftMask = True
ut.ep.ImagingFields.T_tolerance =2.0
ut.ep.ImagingFields.StartReference =1
ut.ep.ImagingFields.NumberOfAverageForReference=4
```

We print the original parameters for thermometry processing

```
[170]: for k in dir(ut.ep.ImagingFields):
        if '_' not in k:
            print(k,getattr(ut.ep.ImagingFields,k))
```

```
Alpha 9.4e-09
AmplitudeCriteriaForRestMotion 25.0
Beta 1.5
CalculateLargeHistory True
CircleSizeFORSNRCoronal 45.0
CorrectionOrder 0
DiscardedPointsInPredictor 100
FrequencyCut 0.8
Gamma 42580000.0
MaxSizeSNRRegion 200.0
NumberNavMessagesToWait 0
NumberOfAverageForReference 4
NumberPointsInterpolateInitialLUT 100
NumberSlicesCoronal 1
OrderForPredictor 5
ROIs 1 C 3 -13 -18
RectSizeFORSNRTransverse 110.0
StartReference 1
TBaseLine 37
TimeBeforeFilterNavigator 10.0
TimeWindowForClassification 11
TimeWindowForFiltering 100
TimeWindowtoKeepInLUT 175.0
UseMotionCompensation False
UseTCoupleForDrift False
UseUserDriftMask True
UserDriftROIs -1 R 70 50 0 -10
```

We process the magntiude and phase data. We also use the CompareTwoOrderedLists to show that the reprocess thermometry is the same as in the original dataset

```
[171]: ut.BatchProcessorFromList(ListDataExample) #Parent class must posses a method
        ↪ directing the processing of entries
```

```
[ ]:
```

We now can plot the different imaging data (magnitude, phase, thermal and supportive mask)

```
[172]: Main=ut.MainObject

def PlotImages(nDynamic, Main,gtitle):
    IMAGES=Main.IMAGES
    plt.figure(figsize=(18,10))
    plt.subplot(2,3,1)
```

```


    p2, p98 = np.
    →percentile(IMAGES['Coronal']['Magnitude'][nDynamic][0]['data'], (2, 98))
    img_rescale = exposure.
    →rescale_intensity(IMAGES['Coronal']['Magnitude'][nDynamic][0]['data'],
    →in_range=(p2, p98))
    plt.imshow(img_rescale, cmap=plt.cm.gray)
    plt.title('Magnitude')
    plt.subplot(2,3,2)
    plt.imshow(IMAGES['Coronal']['Phase'][nDynamic][0]['data'], cmap=plt.cm.gray)
    plt.title('Phase')
    plt.subplot(2,3,3)
    plt.
    →imshow(IMAGES['Coronal']['Temperature'][nDynamic][0]['data'], vmin=40, vmax=55, cmap=plt.
    →cm.jet)
    plt.xlim(100,140)
    plt.ylim(140,100)
    plt.colorbar()
    plt.title('Temperature map')

    plt.subplot(2,3,4)
    plt.
    →imshow(IMAGES['Coronal']['Temperature'][nDynamic][0]['SNR_Mask'], cmap=plt.cm.
    →gray)
    plt.title('SNR mask')
    plt.subplot(2,3,5)
    plt.
    →imshow(IMAGES['Coronal']['Temperature'][nDynamic][0]['SNR_ColdMask'], cmap=plt.
    →cm.gray)
    plt.title('"cold" SNR mask used for drift correction')
    plt.subplot(2,3,6)
    #note that the mask for ROI monitoring is constant accross all image and it
    →is stored in the TemporaryData
    plt.imshow(Main.TemporaryData['Coronal'][0]['MaskAverage']*1.0, cmap=plt.cm.
    →gray)
    plt.title('Mask used for monitoring')
    plt.xlim(100,140)
    plt.ylim(140,100)
    plt.colorbar()
    plt.suptitle(gtitle)


PlotImages(1,Main,'Dynamic=1')#
PlotImages(3,Main,'Dynamic=3')#
PlotImages(5,Main,'Dynamic=5')#


```

Dynamic=1

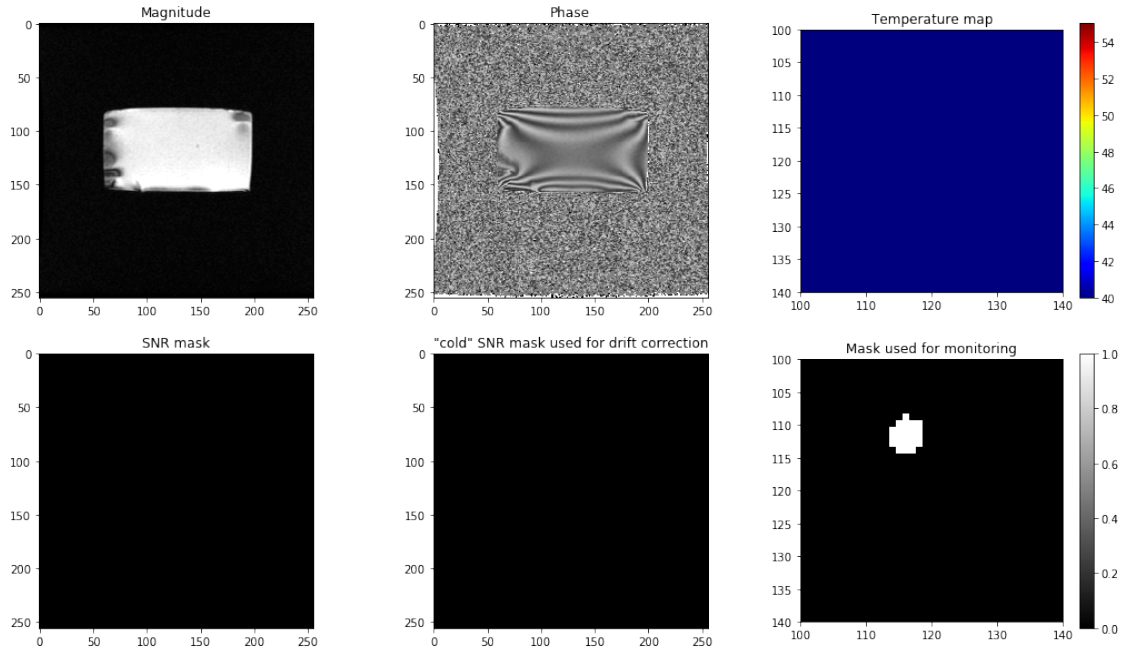

Dynamic=3

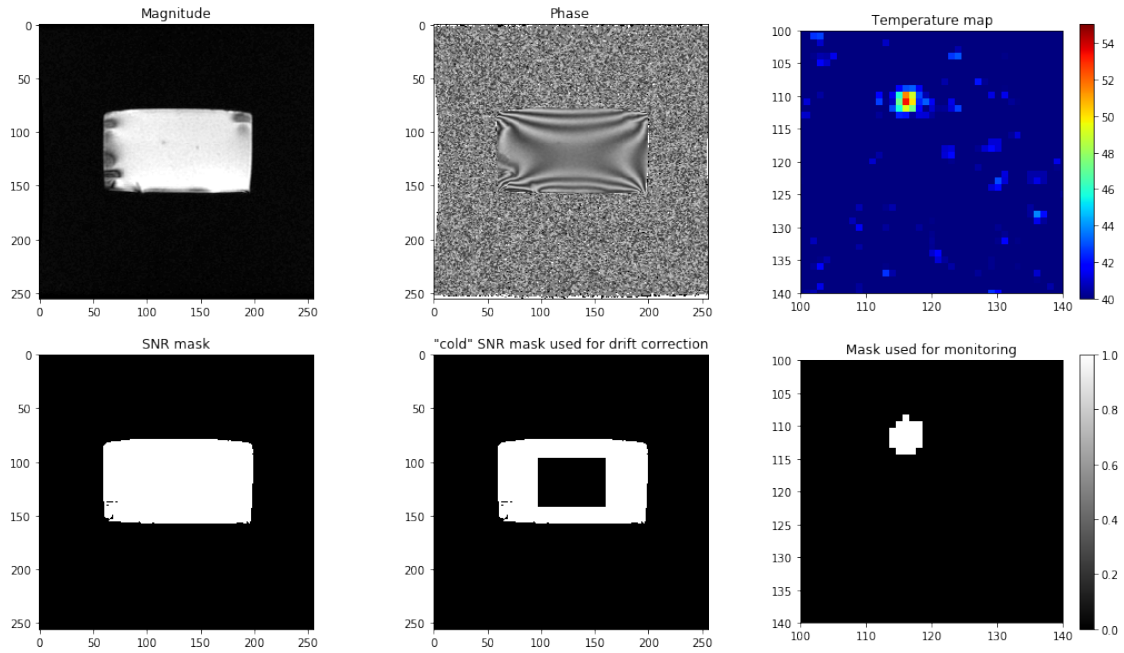

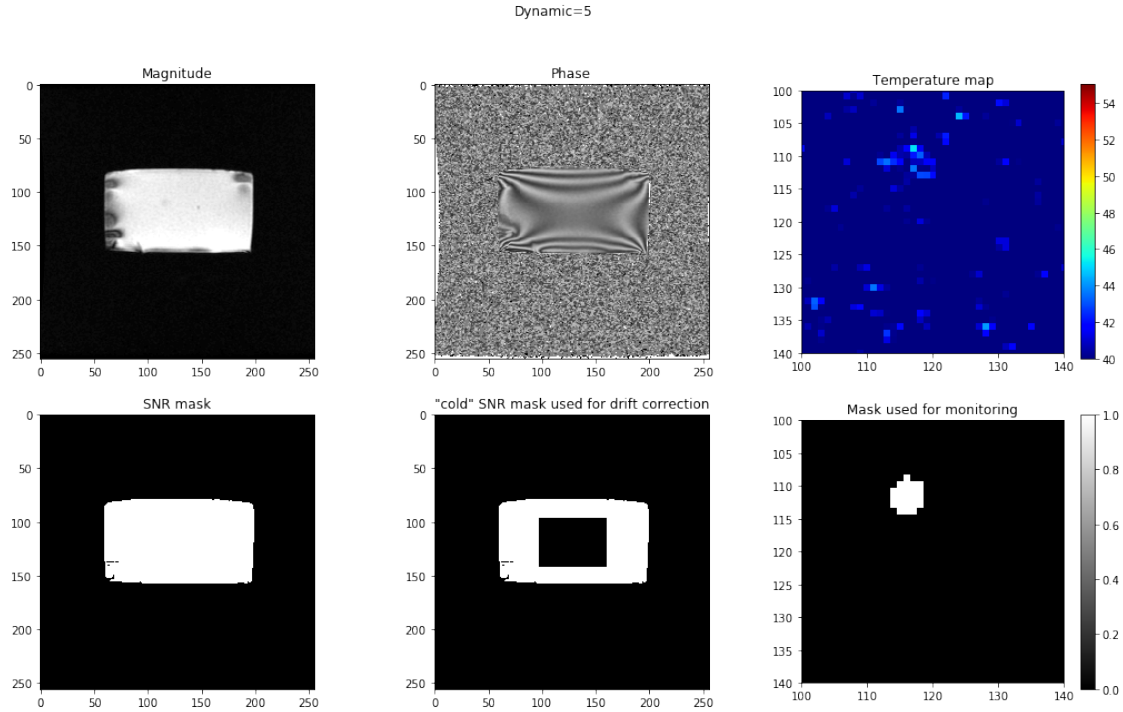

Main.TemporaryData has also the temperature profile over time resulting from the thermometry in the user ROI

```
[173]: def PlotTemporalData(Main):
    timeD=np.array(Main.TemporaryData['Coronal'][0]['TimeTemperature'])
    AvgTemp=np.array(Main.TemporaryData['Coronal'][0]['AvgTemperature'])
    T10=np.array(Main.TemporaryData['Coronal'][0]['T10'])
    T90=np.array(Main.TemporaryData['Coronal'][0]['T90'])
    plt.figure(figsize=(12,6))
    plt.plot(timeD,AvgTemp)
    plt.plot(timeD,T10)
    plt.plot(timeD,T90)
    plt.legend(['Avg. Temperature','T10','T90'])
    plt.xlabel('Time (s)')
    plt.ylabel('Temperature (43$^{\circ}$C)')
```

PlotTemporalData(Main)

